

# thermaeye<sup>®</sup> PLUS

## New System for DRY EYE TREATMENT

Maximum comfort for the patient,  
Maximum safety for the doctor.





**ONE HAND PIECE**

**3 FILTERS**

More Treatments  
All Skin types

The **different wavelengths** available (520, 650 and 800 nm)\* in combination with the **Fractional Square Pulsed Technology**, ensure greater absorption of light by melanin, allowing to work even in the darkest phototypes. Each shot of Thermaeye Plus can be divided up to 3 pulses. The Fractional Square Pulses allow the use higher fluences while protecting the skin.

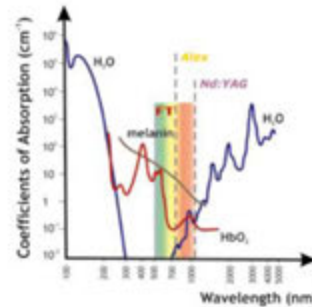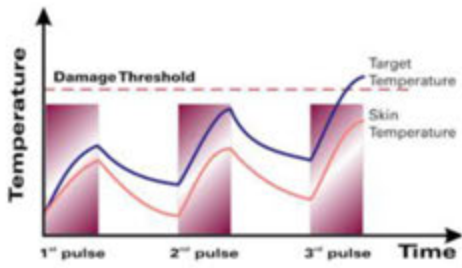

The **new contact cooling system** ensures painless and minimally invasive treatments. This innovative MDS system allows the doctor to eliminate side effects such as redness or burns. **Decreased Risk of Thermal Damage.**

**Very compact and lightweight device.**  
**Easy maintenance.**



Source	Xe Lamp
Spectrum of Emission	520 – 1200 nm (520) 650 – 1200 nm (650) 800 – 1200 nm (SA)
Fluence	Up to 25 J/cm <sup>2</sup>
Pulse Width	From 3 to 8 ms
Delay Between Pulses	From 5 to 50 ms
Number of Pulses	From 1 to 3
Repetition Rate	2 s (min)
Treatment Area	48 mm x 13 mm (6,2 cm <sup>2</sup> )
Integrated Skin Cooling	Contact Peltier System (from 8 to 18°C)
Emission Control	Footswitch and/or Fingerswitch
Electrical Requirements	230Vac or 110Vac 10 A (máx) / 50-60 Hz
Dimensions & Weight	22 cm (H), 45 cm (W), 37 cm (D), 16

Distributed by:



MDS Medical Services S.R.L  
Viale Oberdan, 25  
31100 Treviso (Italy)  
info@thermaeye.com

